


TRAP EASY

FDA
PENDING

AGO PER BIOPSIA OSTEO-MIDOLLARE CON
DISPOSITIVO DI INTRAPPOLAMENTO INTEGRATO

*BONE MARROW BIOPSY NEEDLE WITH
INTEGRATED TRAPPING DEVICE*

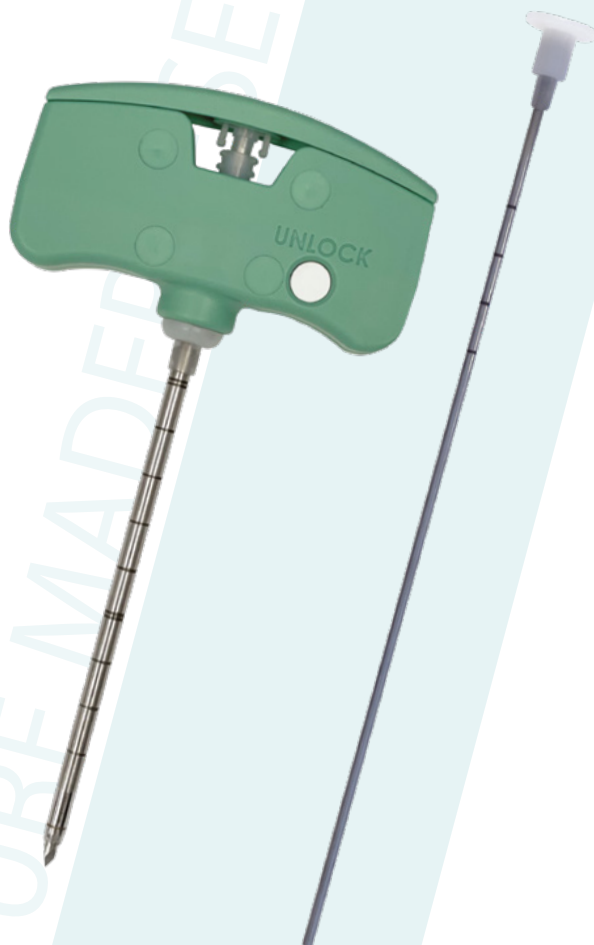

HS[®]
FUTURE MADEPRESENT

TRAP EASY

AGO PER BIOPSIA OSTEO-MIDOLLARE

BONE MARROW BIOPSY NEEDLE

TRAP EASY è un ago per biopsia osseo-midollare che intrappola in modo facile e sicuro il campione e ne assicura il prelievo al 100%.

Caratteristiche del dispositivo:

- facile sistema di attivazione della trappola interna
- campione più grande rispetto ai tradizionali aghi Jamshidi a parità di diametro esterno
- ago cilindrico, centimetrato, con punta a corona per una facile penetrazione del tessuto osseo, con impugnatura ergonomica con connettore luer lock che consente di distribuire la forza su tutto il palmo della mano e non localmente
- mandrino estraibile, perfettamente accoppiato all'ago, con punta piramidale a tre facce
- dispositivo per l'intrappolamento del campione integrato, caratterizzato da un particolare disegno brevettato che consente di imprigionare il frustolo, eliminando la manovra di lussazione
- filo di estrazione a punta smussa, centimetrato, indicato per la valutazione dell'entità del prelievo prima della biopsia e per il recupero del campione
- dispositivo guida di sicurezza indicato per ricoprire la punta dell'ago ed evitare pericoli per l'operatore
- tappo luer lock indicato per chiudere l'ago evitando la fuoriuscita di sangue

TRAP EASY is the a bone-marrow biopsy needle which traps the sample ensuring its complete removal.

TRAP EASY is an bone marrow biopsy needle that easily and safely traps the sample and ensures a 100% sample collection.

Device Features:

- easy internal trap activation system
- larger sample than conventional Jamshidi needles with the same outer diameter
- cylindrical needle with crown tip to easily penetrate the bone, depth marks and ergonomic handle with luer lock connector (it ensures safe grip and manoeuvrability and force distribution on the palm of the hand)
- removable stylet with pyramidal tip
- patented integrated device designed to trap the sample into the needle avoiding the deflection technique
- extraction wire with blunt tip intended for sample size evaluation and sample retrieval
- safety guide intended for needle tip covering to avoid damages for personnel
- luer lock stopper.

CODE	DIAMETER	LENGTH	BOX
TRAPEASY 0810	8G	100mm	10
TRAPEASY 0815	8G	150mm	10
TRAPEASY 1110	11 G	100mm	10
TRAPEASY 1115	11 G	150mm	10